

Handheld OAE-Tympanometry Combination System MAICO ERO•SCAN™ Pro



- OAE and Tympanometry are completed in as little as 20 seconds
- 226 Hz and 1000 Hz Tymp probe tone
- Upgradable to advanced diagnostic capabilities with key codes
- Customizable settings and protocols
- Print to thermal paper, PDF or PC printer
- Optimal patient management database

Handheld OAE-Tympanometry Combination System

MAICO ERO•SCAN™ Pro

DPOAE-TEOAE-Tymp

The MAICO ERO•SCAN™ Pro has been designed for a fast automatic testing of newborn babies, children and adults. It is ideal for ENT-offices, clinics and medical offices.

OAE Tests

All tests begin with AutoStart. The results for up to 10 frequencies are displayed as SNR or Value graph views. The possible outcomes are Pass, Refer, Noisy and Seal.

Tympanometry

The ERO•SCAN™ Pro performs automatic impedance tests with a 226 and 1000 Hz probe tone. The high probe tone is ideal for testing newborns and young children. With the pass criteria "box", results are easy to interpret. Detailed measures are shown on the screen.

Printing

Quick-print to PDF or default PC printer via included software. Results can also be printed by using the thermal printer.

Easy handling


The MAICO ERO•SCAN™ Pro is easy to use, with uncomplicated operation and test analysis. Series Protocols allow multiple tests to be performed by pressing one button.

Patient Management

The optional patient management software/database allows to upload patient names and download test results.

Customizable to individual needs

The models ERO•SCAN™ Pro Standard and Advanced Diagnostic are customizable to individual needs. The front panel programming allows easy change between the basic and a custom test setting. An advanced setup menu allows individual programming of the test parameters.

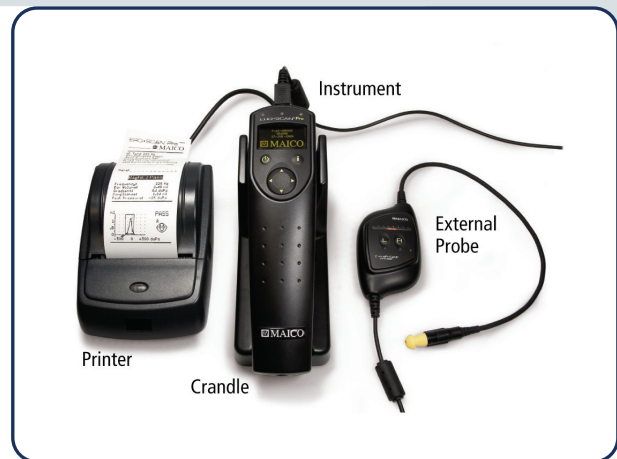

ERO•SCAN™ Pro system with optional external probe

	ERO SCAN Pro DP Standard	ERO SCAN Pro TE Standard	ERO SCAN Pro Advanced Diagnostic
PRODUCT FEATURES			
Internal OAE Probe	Y	Y	Y
External OAE•Probe	Option	Option	Option
External Tymp OAE Probe	Option	Option	Option
DPOAE/TEOAE	DP	TE	DP and TE
SYSTEM FEATURES			
Sequenced Protocols	Y	Y	Y
Quick-print to PDF via PC	Y	Y	Y
Quick-print to PC Printer via PC	Y	Y	Y
Completely Isolated Printer	Y	Y	Y
Completely portable	Y	Y	Y
Memory (no. tests) Maximum	350	350	350
Save tests by patient no.	Y	Y	Y
Prints graphical and numeric data	Y	Y	Y
Date & Time on printout	Y	Y	Y
Tests all ages	Y	Y	Y
Patient Management Database	Option	Option	Option
OAE FEATURES			
OAE Frequency Range (Hz)	1.5-12DP	0.7-4 TE	1.5-12DP 0.7-4 TE
Maximum no. of Test Frequencies or Bands Reported	10 DP	10 TE	10 DP 10 TE
Default Pass/Refer	Y	Y	Y
No. of OAE Test Protocols	20	20	20 DP 20 TE
Completely Isolated Printer "optional" (does not require external line source)	Y	Y	Y
Customizable Pass Criteria	Y	Y	Y
Customizable Test Protocols	Y	Y	Y
Customizable by frequency/range: stimulus intensity, P1/P2 ratio, averaging time, SNR for pass, minimum OAE amplitude	Y	Y	Y
Selectable Averaging Time	Y	Y	Y
TYMPANOMETRY OPTION			
Probe tone 226 Hz	Y	Y	Y
Probe tone 1000 Hz	Y	Y	Y
Default PASS/REFER indicator & box	Y	Y	Y
Show/Hide PASS/REFER indicator & box	Y	Y	Y
Customizable Pass Criteria	Y	Y	Y
No. of Tymp Protocols	10	10	10

ERO•SCAN™ is a trademark of ETYMOTIC RESEARCH



MAICO Diagnostic GmbH
Salzufer 13/14, 10587 Berlin / Germany
Tel.: + 49 30/70 71 46 50
Fax: + 49 30/70 71 46 99

e-mail: export@maico.biz
Internet: www.maico.biz